sinus

Body

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Contents English

Operating instructions	18
1. Introduction	19
2. Description of the device	19
3. Contraindications	20
4. Warranty	21
5. Operating and handling	22
6. Technical specifications	26

Operating instructions sinus Body

Before putting the device on for the first time, you need to be aware the restrictions for wearing the sinus Body listed under section 3!

The living organism has been evolving and adapting to Mother Earth for millennia. A healthy person is in resonance with the Earth.

However, in recent decades technical developments are overtaking us at a rapid pace. The speed of technological evolution is so high that the human organism doesn't have enough time to adapt to the respective innovations and pressures. Moreover, there are patterns of behaviour that are not particularly conducive to the body's own capacity to regenerate (staying in buildings shielded from geomagnetism, too little time in nature, longer periods in Faraday cages [cars, public transport, concrete], improper diet, stress, etc. In addition, the natural geomagnetic field continues to decline (the earth's magnetic field has lost about 50% of its power in the past 500 years). This can cause physical and mental stress, which means we lose a lot of energy. The body needs this energy to maintain our vitality. This is where the sinus Body comes in.

The sinus Body is a pulsed magnetic field that supplies the body's own cells with the frequency of the earth of 8 hertz (Schumann frequency, 1952). The pulsed field enhances a person's own field. As a result, the body is able to better process stress (bioresilience). The body's cells are given the opportunity to regenerate on a natural level, thus increasing vitality.

Supplemented by an innovative form of information field exposure (based on the principles of bioresonance regulation), the so-called electrosmog of the sinus Body itself becomes aligned

and no longer creates interference.

1. Introduction

The sinus Body is a technical device for generating magnetic field pulses that match those of nature.

Please read this manual carefully for instructions on handling the device. Using this device is not a replacement for a healthy lifestyle, which includes sufficient sleep, outdoor exercise, and a healthy diet. As a trademarked wellness device, the sinus Body is not a replacement for treatment by a responsible doctor or therapist.

Further information and literature references for a better understanding of the effects on the human organism and the effects of pulsed magnetic field technology can be obtained from your consultant on request.

The normal scope of delivery includes:

- The sinus Body magnetic field device
- Magnet for activating the device
- Lanyard and jewellery box
- Induction charger
- Operating manual

2. Description of the device

The sinus Body is a miniaturised, pulsed magnetic field device. It has battery-powered electronics that produce a pulsed magnetic field matching that of nature, at a frequency of 8 Hz (hertz). This means it induces 8 pulses per second. This is the magnetic field that should be transmitted to humans through natural geomagnetism.

The sinus Body is also informed by a specially developed field exposure technology for detecting interference frequencies and for neutralising them via destructive interference. This works on the principle of bioresonance regulation.

In addition to the magnetic impulses and the counter frequencies, the sinus Body also emits coloured light impulses. The sinus Body is worn on the chest with the light turned towards the body. Thus, the organism is stimulated from here via the main meridian (Ren-Mai, the conception vessel).

3. When the sinus Body should not be worn

People with the following limitations should not use pulsed magnetic field systems, including the sinus Body (or only after consultation with a physician):

- Those who have cardiac pacemakers and other electronic implants (e.g. insulin pumps).
- People with severe angina pectoris.
- People with acute bleeding.
- People with severe fungal disorders.
- People who suffer from epilepsy.
- People with severe cardiac arrhythmias (grades 3 and 4).
- People who suffer from acute tuberculosis.
- People with hyperthyroidism.
- People with a fever (over 39° C).

During pregnancy, magnetic therapy should only be applied in exceptions and after professional consultation.

4. Warranty

The warranty period for a device is 24 months from the date of delivery. Wear parts, consumables, case surfaces, necklaces and eyelets are generally excluded from warranty claims.


A warranty claim expires:

- in the event of damage due to improper operation or inappropriate use of the sinus Body;
- in the event of damage due to mechanical influences caused by the customer, such as falling, crushing, contact with water, etc.;
- in the event of tampering with the device's technology or use of wear, consumable or spare parts which are not original parts or have not been declared suitable as described in this manual;
- if changes, extensions or repairs are made to the device by unauthorised people;
- if this operating manual is not complied with when using the device;
- if seal stamps, lead seals or corresponding securing fasteners have been loosened, removed or separated.

The start of the warranty period is based on the official delivery date of the company i-like Metaphysik and the associated delivery to the purchaser. The official delivery note or the official invoice issued by the company i-like Metaphysik is required in order to assert the warranty.

5. Operating and handling

Settings for getting started

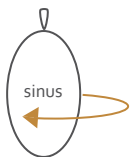1. The sinus Body is delivered turned off and charged. Take the sinus Body out of the jewellery box and place it on the table. Turn the sinus Body onto its back so that the light is pointing up.

2. To turn on the sinus Body, take the supplied magnet and hold it to the lower left end of the device housing. Look at the small light in the case (above, the small hole near the eyelet). When it activates, you will see a light flashing every 10-12 seconds. The sinus Body is now turned on.


3. Now thread the supplied lanyard (or your own necklace*) through the oval eyelet on the housing. The sinus Body is worn with a necklace around your neck, with the light turned to your body, directly on the skin. Please note, special care is required if you pull the closed necklace over your head. Breaking of the necklace or the housing eyelet is not covered by the warranty (material defect liability).

*You can of course replace the lanyard with another necklace. If you use a metallic necklace, please make sure that it is galvanically interrupted once. This means that you interrupt the metallic necklace laced through the device eyelet with a non-metallic chain link (e.g. in the vicinity of the closure), or you can use a completely non-metallic necklace, like one made of rubber, leather, silk cloth, or textile.

Important information

You should drink enough water (non-carbonated) or tea per day, so that the organism can respond as well as possible to your optimised self field. You can also wear the sinus Body at intervals (e.g. during the day only by hourly intervals) to get used to it. Most wearers report that they prefer taking the sinus Body off at night.

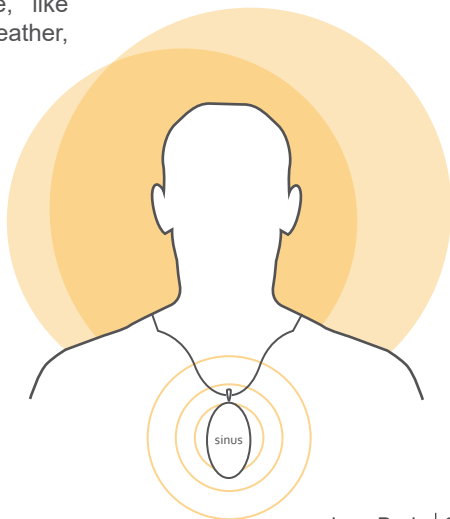Turning it on

The device is off if the light is not flashing. To turn on the sinus Body, turn it so that the light is visible. Take the supplied magnet and hold it to the lower left end of the device housing until a light appears on the LED. The light will then flash in a cycle of approximately 10-12 seconds.

Turning it off

The device is on if the light is flashing in an approximately 10 to 12-second cycle. To turn the sinus Body off, turn it so that the light is visible. Take the supplied magnet and hold it to the lower left end of the device housing for at least 3 seconds, until the light flashes red several times. The light will no longer be visible after that. The sinus Body is now turned off.

Charging the battery

If the light is dim or is not lighting up even if it is switched on, the battery is running low or the device has been switched off.

Please recharge the battery using an induction charger. Place the sinus Body on the charger so that the light is pointing upwards. The charging time depends on the charger. The charger will emit a signal when the battery is fully charged (the sinus Body itself will not alert you that it is charged). Use the supplied original charger or an induction charger that is commonly used for mobile phones.

If the sinus Body cannot be charged properly, this is almost always due to the charger or the power source it is plugged into (AC adapter, cell phone charger, etc.). The operating time per battery charge is about 2-3 weeks. If the device switches off due to insufficient battery charge, you need to restart the sinus Body after recharging!

Note: magnetic field impulses are only produced for as long as the power electronics are supplied with sufficient voltage from the battery. We therefore recommend charging the device once per week.

Note on the device's water resistance

The sinus Body's housing is not waterproof, so you must not shower, bathe or do water sports while wearing the device.

Cleaning

Clean the outside of the device housing with a soft, clean cloth. Please do not use „aggressive“ chemical household cleaners.

Disposal

We would ask you not to dispose of the sinus Body in residual waste, but to instead bring the device to a suitable recycling centre or collection site. The device is registered with WEEE – its disposal is paid (recyclable electronics).

Note regarding scientific acceptance

As a manufacturer, we comply with the obligation to make no healing claims with respect to the sinus Body. The statements made and information on the effects of the device and the pulsed magnetic field are derived from source-based scientific studies and user reports.

We dissociate ourselves from any healing claims and promises of a cure, as well as from claims of medical effectiveness in the sense of the Therapeutic Products Act (TPA) and the Unfair Competition Act (UCA).

6. Technical specifications

Material housing:	birch wood
Housing construction:	three-dimensional, elliptical, oblong
Housing dimensions (approx.):	elliptical, (W/D/H) 30 mm x 20 mm x 60 mm
Device weight:	approximately 15 g
Operating temperature:	10°C - 40°C
Storage temperature:	0°C - 70°C
Contents:	electronic pulse generator
Type of coupling:	inductive coupling of a patented magnetic emitter
Field strength:	pulse value 2000 μ T
Frequency:	pulsed frequency pattern as a sinus
Supply voltage:	battery
Battery type:	LiPo401230
Rated capacity:	105 mAh
Rated voltage:	3,7 V
Charging terminal voltage:	4,2 V

Discharge voltage:	3,0 V
Charging current:	75 mAh
Magnetic impulses:	every 10 sec. over 1 sec. duration in the cycle of 8 Hz
Light:	green LED in 10 sec intervals when functioning red and green LED for switching states
Conformity:	CE as a wellness device with declaration of conformity
EMC:	electromagnetic compatibility (EMC) certified
WEEE:	registered
Trademarks:	“sinus i-like” as a registered trademark by i-like Metaphysik GmbH, CH-9445 Rebstein
Country of origin:	Switzerland + Germany

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