sinus Body

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Contents English

User manual

1. Introduction	18
2. Description of the device	18
3. Contraindications	19
4. Important safety instructions	20
5. Warranty	21
6. Operating and handling	22
7. Technical specifications	26

User manual sinus Body

1. Introduction

Before wearing the device for the first time, you need to be aware of the restrictions for wearing the sinus Body listed under section 3!

The sinus Body is a technical device for generating magnetic field pulses in the form of a sinus wave that match those of nature. Please read this manual carefully for instructions on handling the device. Using this device is not a replacement for a healthy lifestyle, which includes sufficient sleep, outdoor exercise, and a healthy diet. As a trademarked wellness device, the sinus Body is not a replacement for treatment by a responsible doctor or therapist.

Further information and literature references for a better understanding of the effects on the human organism and the effects of pulsed magnetic field technology can be obtained from your consultant on request.

The normal scope of delivery includes:

- The sinus Body magnetic field device
- Lanyard and jewellery box
- USB charging cable
- USB-C cover
- Operating manual

2. Description of the device

The sinus Body is a miniaturised, pulsed geomagnetic field device. It has battery-powered electronics that produce a pulsed magnetic field matching that of nature, at a frequency of 8 Hz (hertz). This means it generates 8 pulses per second. This is the magnetic field that should be transmitted to humans through natural geomagnetism.

The sinus Body is also informed by a specially developed field exposure technology to make vital frequencies available to the body. In this way, the body's resilience can be strengthened. This works on the principle of bioresistance regulation.

In addition to the sinus frequency, the sinus Body also emits coloured light impulses. The sinus Body is worn on the chest (at the heart chakra) with the light turned towards the body.

3. When should the sinus Body not be worn?

People with the following limitations should not use pulsed magnetic field systems, including the sinus Body (or only after consultation with a physician):

- those who have cardiac pacemakers and other electronic implants (e.g. insulin pumps)
- People with severe angina pectoris
- People with acute bleeding
- People with severe fungal disorders
- People who suffer from epilepsy
- People with severe cardiac arrhythmias (grades 3 and 4)
- People who suffer from acute tuberculosis
- People with hyperthyroidism
- People with a fever (over 39 °C)

During pregnancy, use of magnetic field systems should only be applied in exceptions and after professional consultation.

4. Important safety instructions

With this device, you acquire a modern device that contains a Li-Ion battery. The lifespan of the device essentially depends on the condition and handling of the built-in battery.

If the green LED starts to flash twice every 12 seconds, the device will soon need to be recharged. If the device is not recharged in time, it will be deactivated automatically.

The device should not be stored for more than 30 days in an empty state without recharging. Otherwise, you risk damaging the battery.

Make sure that the charging cable is plugged into the device correctly and without massive pressure. When plugging in, ensure that the plug is guided straight. Do not tilt it! Only use the charging cable supplied.

Please always recharge the device after use. The device should not be permanently connected to the charger as this accelerates the ageing of the battery.

Do not expose to any heat source.

Do not store on top of a heat source.

Do not throw into fire.

Do not store in your car.

Please note the storage instructions.

Do not dispose of with household waste.

5. Warranty

The warranty period for a device is 24 months from the date of delivery. Wear parts, consumables, case surfaces, necklaces and eyelets are generally excluded from warranty claims. A warranty claim expires:

- in the event of damage due to improper operation or inappropriate use of the sinus Body;
- in the event of damage due to mechanical influences caused by the customer, such as falling, crushing, contact with water, etc.;
- damage, failure, loss or deterioration due to ordinary wear and tear, from abuse, misuse, negligence, accident, faulty installation, placement and operation, improper maintenance, unauthorized repair, exposure to unreasonable temperatures or conditions;
- in case of incorrect installation or loss of accessories;
- in the event of tampering with the device's technology or use of wear parts, consumables or spare parts which are not original or have not been declared suitable by this manual;
- if changes, extensions or repairs are made by unauthorised people;
- if this operating manual is not complied with when using the device;
- if stamps, lead seals or corresponding securing fasteners have been loosened, removed or separated;
- damage resulting from inadequate packaging or carrier mishandling of returned products.


The start of the warranty period is based on the official delivery date of the i-like Metaphysik company and the associated delivery to the purchaser. The official delivery note or the official invoice issued by the i-like Metaphysik company is required in order to assert the warranty.

6. Operating and handling

Settings for getting started

1. The sinus Body must be charged before it is worn for the first time. Take the supplied charging cable and plug it into the sinus Body. Then plug the other end of the cable into a suitable USB socket of a mains plug or device with a power supply (e.g. computer).

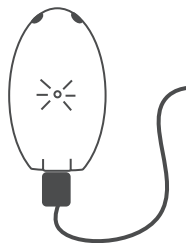2. As soon as the charging function starts, the light on the back of the housing flashes red once per second. The sinus Body is now in charging mode. When the battery is charged, the light signal changes to green and flashes every 12 seconds. The sinus Body is now charged and switched on.

3. Now thread the supplied lanyard (or your own necklace*) through the oval eyelet on the housing. The sinus Body is worn with a necklace around your neck, with the light turned to your body, directly on the skin. Please note, special care is required if you pull the closed necklace over your head. Breaking of the necklace or the housing eyelet is not covered by the warranty.

4. Keep connections dry and clean. To protect the USB socket of the sinus device from dust, dirt and moisture, please use the supplied cover for the USB socket. Simply insert the cover straight into the USB socket and only remove it when the device is to be charged.

*You can of course replace the lanyard with another necklace. If you use a metallic necklace, please make sure that it is galvanically interrupted once. This means that you interrupt the metallic necklace laced through the device's eyelet with a non-metallic chain link (e.g. in the vicinity of the closure), or you can use a completely non-metallic necklace, like one made of rubber, leather, silk, or textile.

Important information

You should drink enough water (non-carbonated) or tea per day, so that the organism can respond as well as possible to your optimised self field. You can also wear the sinus Body at intervals (e.g. during the day only by hourly intervals) to get used to it. Most wearers report that they prefer taking the sinus Body off at night.

Switching on

The device is switched off when no light signal is visible. In this case the battery is empty (deep discharge protection for the battery is integrated). To switch on the sinus Body, plug the supplied charging cable into the sinus Body. Then plug the other end of the cable into a suitable USB socket of a mains plug or device with a power supply (e.g. computer). As soon as the sinus Body charging function starts, the light on the back of the housing will flash red once per second. When the battery is charged, the light signal changes to green and flashes every 12 seconds. The sinus Body is now charged and switched on.

Switching off

The device is switched on when a light signal lights up approx. every 12 seconds. The sinus Body remains switched on until the battery is used up. Then the sinus Body switches off automatically.

Battery charge and light signals

When the sinus Body is switched on and charged, a green light signal lights up **once** every 12 seconds on the back of the housing. If the battery level is low, the light signal flashes **twice** every 12 seconds. If the sinus Body is not charged in spite of a low battery level, it switches off automatically when the battery is empty (integrated deep discharge protection for the battery).

Please charge the battery with the charging cable supplied. To do this, plug one end of the cable into the sinus Body and the other end into a suitable USB socket of a mains plug or device with supply current (e.g. computer). As soon as the charging function starts, the light on the back of the housing flashes red once per second. The sinus body is now in charging mode. When the battery is charged, the light signal changes to green and flashes every 12 seconds. The sinus body is now charged and switched on.

Note: Magnetic field pulses are only produced as long as the electronics can handle the voltage of the battery. We therefore recommend charging the device every 1–2 weeks.

Note on the water resistance of the device

The sinus Body's housing is not waterproof, so you must not shower, bathe or do water sports while wearing the device.

Cleaning

The housing must be cleaned regularly. Clean the outside of the device's housing with a soft, clean cloth. Please do not use harsh chemical household cleaners.

Disposal

We would ask you not to dispose of the sinus Body in general waste, but to instead bring the device to a suitable recycling centre or collection site.

Note regarding scientific acceptance

As a manufacturer, we comply with the obligation to make no healing claims with respect to the sinus Body. The statements made and information on the effects of the device and the pulsed magnetic field are derived from source-based scientific studies and user reports.

We dissociate ourselves from any healing claims and promises of a cure, as well as from claims of medical effectiveness in the sense of the Therapeutic Products Act (TPA) and the Unfair Competition Act (UCA).

7. Technical specifications

Housing material:	TBD
Housing construction:	Three-dimensional, elliptical, elongated
Housing dimensions (approx.):	Elliptical (W/D/H) 27 mm x 16 mm x 46 mm
Device weight:	Approx. 15g
Operating temperature:	10 °C – 40 °C
Storage temperature:	0 °C – 70 °C
Content:	Electronic pulse generator
Frequency:	Pulsating frequency pattern as a sine wave
Supply voltage:	Battery
Battery Type:	Li-Ion
Nominal capacity:	40 mAh
Nominated voltage:	3.7 V
Power consumption charging mode at USB-C socket:	max. 40 mA
Magnetic impulses:	Continuous pulse at a rate of 8 Hz

Light signals:	LED green in 12 sec. flashing indicates device is on. Red and green LEDs for switching states.
Conformity:	60335-1 el. Safety CE Wellness device
EMV:	2014/30/EU EMC directive
Country of origin:	Switzerland
Trademarks:	“sinus i-like” as a registered trademark of: i-like Metaphysik GmbH CH-9445 Rebstein

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